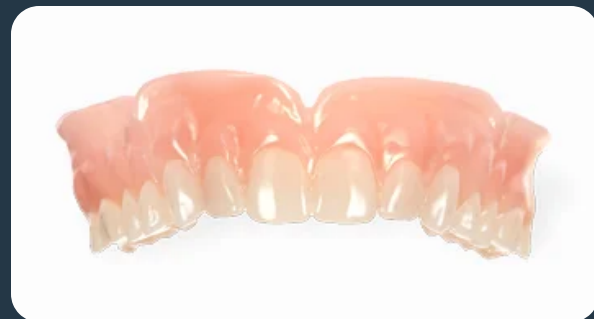


FULL & PARTIAL CASES

DIGITAL DENTURES

Fewer appointments. Less chairside time. A superior fit for your edentulous patients — designed digitally by the same technicians who answer your phone.

3D-printed in our Bartlett, TN lab — precision you can seat.

**2–3****Appointments per case**

vs. up to 5 the conventional way

5 Days**In-lab turnaround**

Made in Bartlett, TN · shipped overnight

5 Years**Your file stays on record**

Remakes & spares from the original scan — no new appointment

THE R-DENT DIFFERENCE

One intraoral scan replaces the extra impression appointments of a conventional workflow — no new technique to learn. Our technicians design the case in the latest CAD software and 3D-print the base and teeth for a denture that fits right the first time. The digital file stays on record for five years, so remakes and spares start from the original record, not a new appointment.

Two Ways to Deliver a Digital Denture

2 APPOINTMENTS**\$149–\$249 / arch**

Reference Denture

For patients with an existing denture. We use their current denture as the reference to capture the impression and bite — no try-in, and they keep their denture. Also the fastest path to an exact copy or spare.

- 1 Scan & order:** scan the existing denture extraorally (a wash impression adds detail) and scan the bite intraorally.
- 2 Seat:** deliver the finished denture. Done.

3 APPOINTMENTS**\$199–\$299 / arch**

New Denture

No existing appliance

For edentulous patients with no current denture, or when esthetics, bite, or VDO must change. We capture the record with a Glenn Fit tray — no bite rims — then a printed try-in lets you and the patient approve the result before the final.

- 1 Scan & order:** seat a **Glenn Fit tray** for the initial impression, scan the arches, and request a printed try-in.
- 2 Register bite / try-in:** scan the bite and approve esthetics and fit.
- 3 Seat:** deliver the finished denture.



Glenn Fit trays keep this at 3 appointments — a custom tray adds one. Scan to order your kit.

Why Dentists Choose R-Dent

- ✓

Real people, not a phone tree. Talk directly to the technicians building your case.
- ✓

Owner-operated since 1978. Family-owned, not a venture-backed platform.
- ✓

Scanner-agnostic. Use the intraoral scanner you already own — no lock-in, no minimum spend.
- ✓

Made in the USA. Fabricated in Bartlett, TN and shipped overnight via UPS.

The R-Dent Difference at a Glance

	R-Dent	Big Digital-Only Labs
Support	Real people who answer the phone — no ticket queue	Chat bots and call-center support tiers
Scanner	Scanner-agnostic — use the one you already own	Often tied to a proprietary scanner and minimum spend
Digital file	On record for 5 years — reprints and spares without a new scan	Retention policies vary by platform
Ownership	Family-owned since 1978, not a venture-backed platform	Scaled nationwide, VC-funded

Using the Glenn Denture Tray

The Glenn Fit tray captures the initial impression for a New Denture case in one appointment — no bite rims, no custom tray. Starter Pack \$159 covers all sizes, both arches; one case recovers the kit cost 2x.



SCAN: BUY TRAYS



SCAN: HOW TO USE

Ready to send a case?

Scan the arch, submit your digital file, and let our technicians handle the rest.

CALL 901.372.8020

TEXT 901.201.5336

rdentlab.com/send-a-case