

## Clinical Technique Spotlight

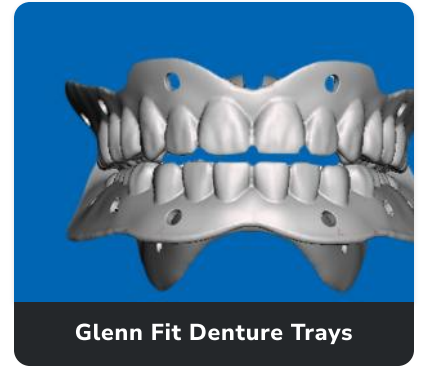
# COMPLETE DENTURES IN 2 APPOINTMENTS

A thermoplastic edentulous impression tray with teeth built in — it combines the master impression, jaw relation record, and an esthetic preview into a single first visit, replacing the baseplate-and-wax-rim appointment entirely.

S / M / L • UPPER &amp; LOWER

ANALOG + DIGITAL COMPATIBLE

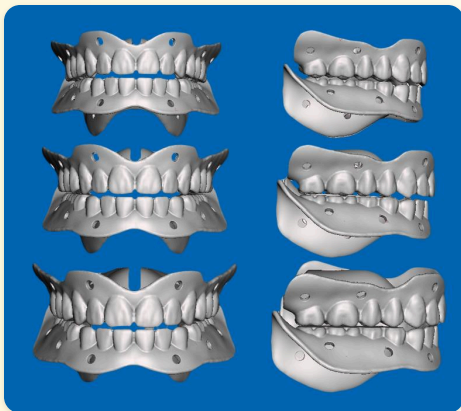
Developed by Dr. N. Cory Glenn, DDS • VP of Technology, Blue Sky Bio

**5 →  
2**

Appointments

**~\$272+**

Net gain per case

**~\$259**Appt slot  
recovered (ADA  
HPI)**\$29**Tray 2-  
pack**THE FULL KIT****Three Sizes, Both Arches**

Small, Medium, and Large trays, front and side profile, each pre-formed with anatomic teeth and retention holes for a secure impression material bond.

**\$29**

Single Arch 2-Pack • 1 size

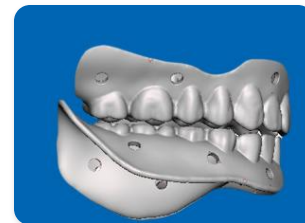
**\$159**

Starter Pack • 12 trays, all sizes



## How It Works — One Visit

Appointment 1 replaces the primary impression, wax-rim try-in, and master impression — all in 30–40 minutes.



1

### Select Tray Size

Measure nose width with calipers — matches canine-to-canine width. Medium fits most adults.

2

### Heat & Pre-Adapt

Submerge in near-boiling water ~15 sec. Widen or narrow arch, bend flanges to fit. Cool to set.

3

### Line & Set VDO

Soften thermoplastic beads, line tray interior. Seat and align to Fox plane and ala-tragus line.

4

### Trim & Border Mold

Trim long borders, extend posterior border to hamular notches. Verify peripheral seal.

5

### PVS Impression

Apply tray adhesive. Heavy body on borders, light body internally. Repeat for opposing arch.

6

### Bite + Photo + Send

Confirm centric relation, take a straight-on photo for tooth reference, and send to the lab.

### Guaranteed ROI — Every Case

**~\$259** appt slot + **\$35** lab fee — **~\$22** tray cost = **~\$272+** net gain

Eliminating the dedicated wax-rim try-in frees a 30–60 minute appointment slot for a billable procedure — per ADA Health Policy Institute data, worth ~\$259 on average. Guaranteed, every single case, regardless of fee level.

## Ready to try the Glenn Denture Tray?

Starter Pack \$159 — all sizes, both arches. One case recovers the kit cost 2x over.

Buy Trays at [blueskybio.com](https://blueskybio.com)



Order Trays  
[blueskybio.com](https://blueskybio.com)



Watch Tutorial  
42-min video

\* ADA HPI Survey of Dental Practice, 2024. Appointment slot value ~\$259 and tray cost ~\$22/arch pair are averages; individual results vary. Glenn Denture Tray developed by Dr. N. Cory Glenn, DDS — VP of Technology, Blue Sky Bio · [coryglenn.org](https://coryglenn.org). Have questions about digital denture cases? R-Dent is here to help.